

Donning Enhanced PPE – Aerosol Generating Medical Procedures (AGMP)

1) PREPARE

- Remove watches, jewelry and any personal items. Ensure ID lanyard is covered or removed
- Secure hair back
- Be clean shaven

2) ASSEMBLE SUPPLIES

- Long sleeved gown, N95 Respirator, Face Shield, extended cuff gloves, bouffant, hand hygiene supplies

3) PERFORM HAND HYGIENE

4) DON GOWN


Secure neck ties



Secure waist ties



Wear gown close to neck

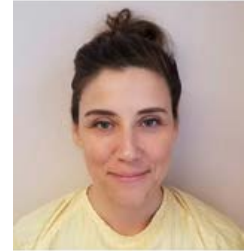

5) DON N95 RESPIRATOR, BOUFFANT & FACE SHIELD

Don N95 Respirator
Straps DO NOT cross



Don Bouffant
Cover all visible hair



Mold N95 Respirator to nose
Fit mask under chin
Test for seal



Add face shield
Adjust to protect lower neck



5) DON GLOVES

Gloves over cuff



READY

Enhanced PPE with N95 Respirator and Face Shield

